

Remote learning

Screen mirroring: connection, mirroring, data acquisition

Connection to additional end devices: 3 options


WLAN connection



Direct LAN connection




Integration in customer's own network

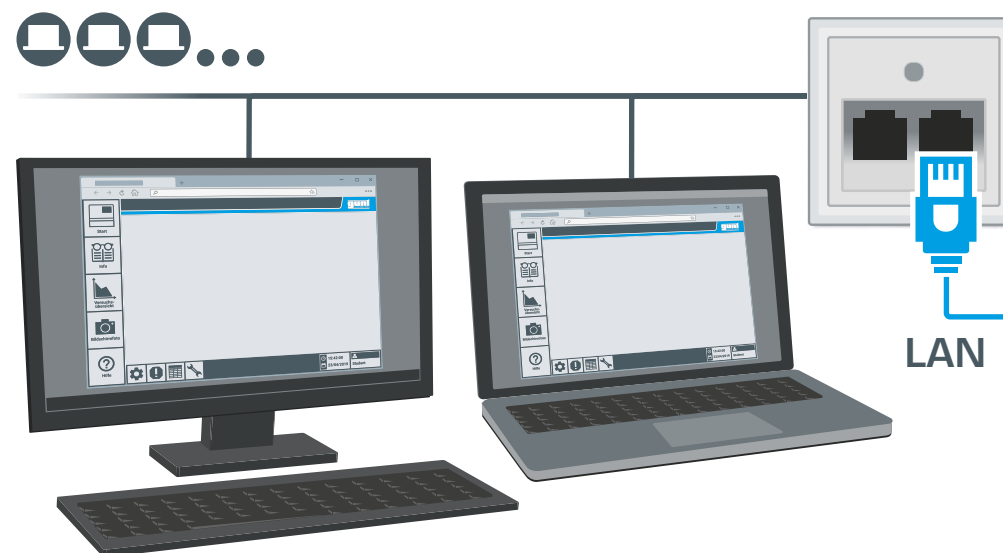

Touch screen with integrated PLC


Internal memory or USB flash drive for measured values

WLAN router integrated in the experimental unit

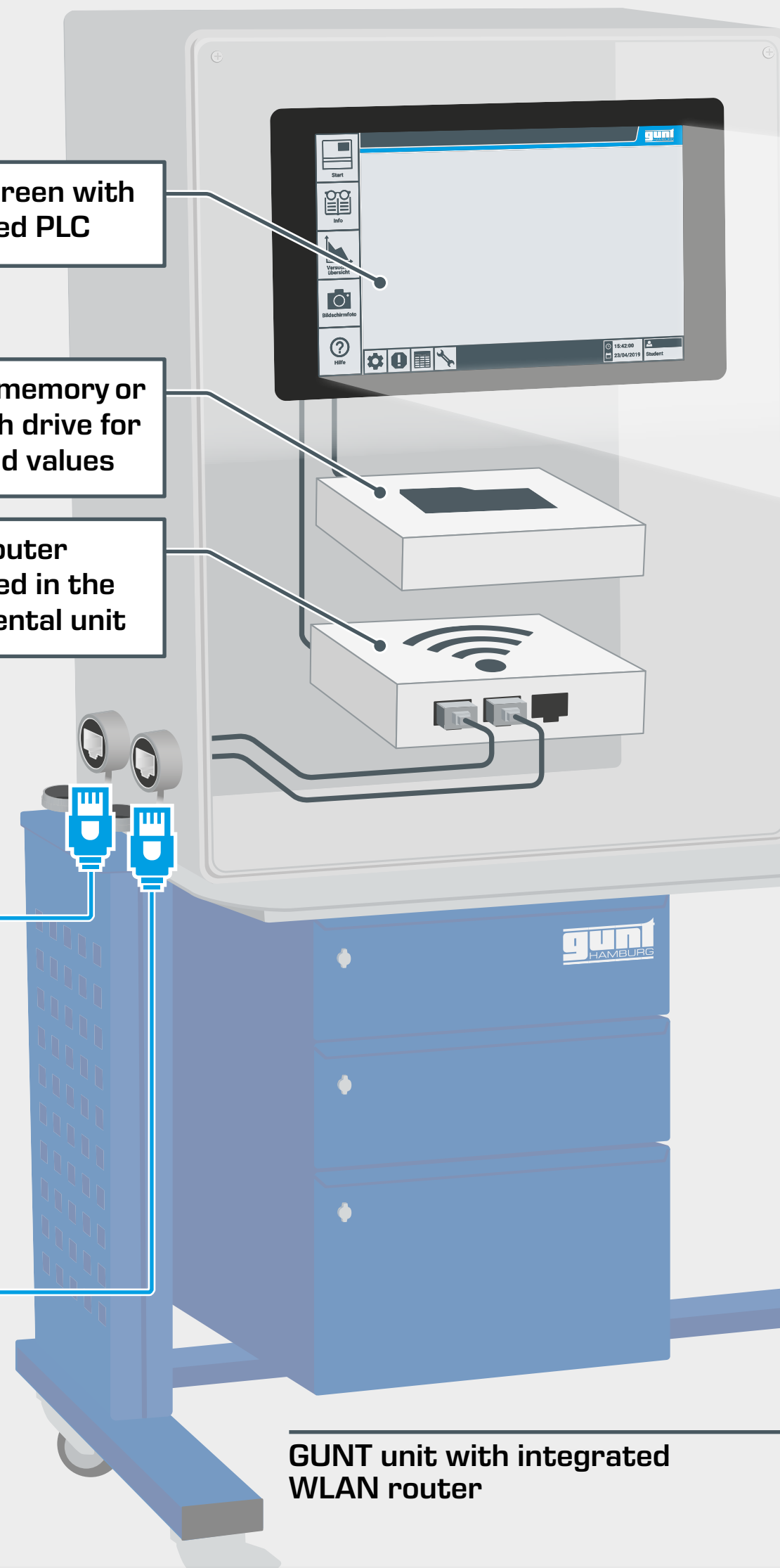

GUNT unit with integrated WLAN router


Mirroring of the user interface

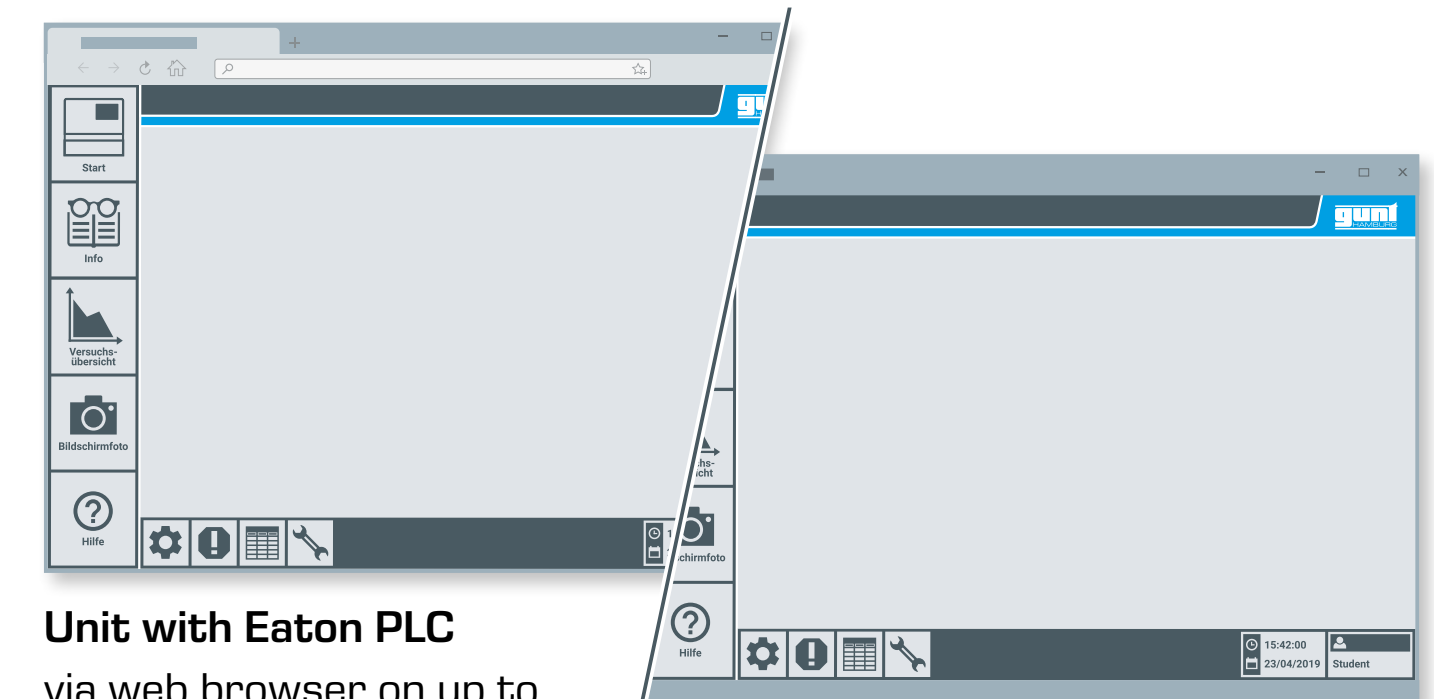

Unit with Eaton PLC
via web browser on up to 5 additional end devices

Unit with Weintek PLC
via installed cMT Viewer on up to 10 additional end devices

Transfer measured values

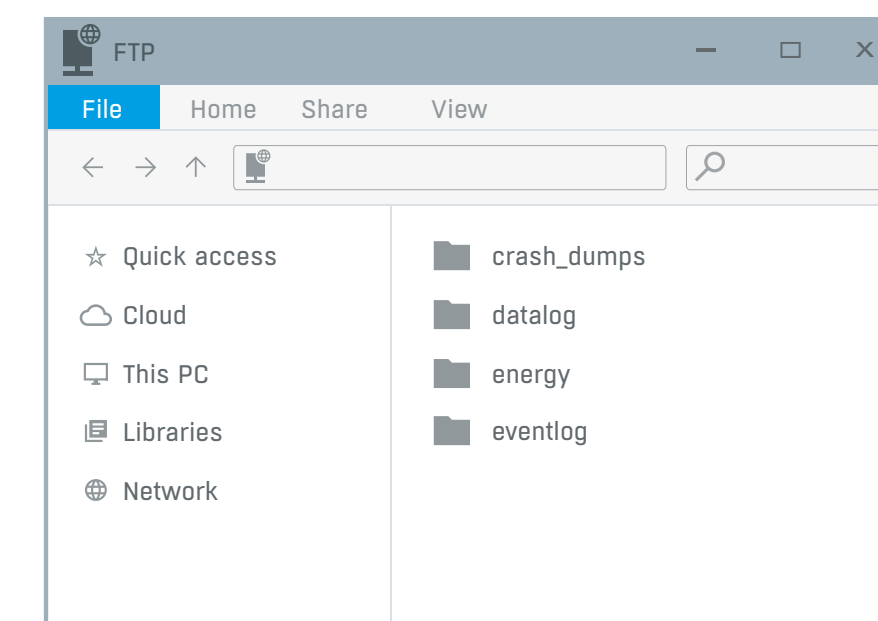

By means of the IP address of the experimental unit, it is possible to access measured values that were stored internally or via USB flash drive during the experiment using an FTP server.